


Medical Procedure Masks

Medical Procedure Masks are constructed of a three-layer, nonwoven fabric featuring a spunbond polypropylene on the inner and outer layers and a meltblown polypropylene middle layer. Pleated Earloops, One Size Fits Most.

Meets the ASTM requirements for Level 1 barrier protection.

Bacterial and particulate filtration efficiencies greater than 95 percent with fluid resistance at 80 mm Hg.

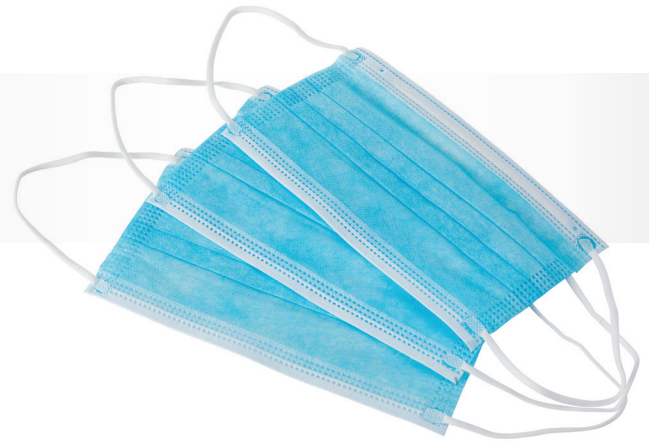

Features & Benefits:

- One Size Fits Most
- ASTM Level 1 Medical Face Masks
- Cool and breathable SMS material with pleated design
- Single-strand, elastic ear loops designed to be low-profile and stretch to enhance comfort and fit with less irritation
- Not made with natural rubber latex (Latex Free)

Product	Reorder #	Units
Medical Procedure Masks	PC622B	50 Masks/Box, 20 Boxes/Case

Visit Us Online at:
www.ProCureProducts.com